

Rivarp

MEDICAL

MEDITRANS Contact Info



Jakarta : Jl Munggang no. 86-87 Balekambang
Kramat Jati, Jakarta Timur, 13530
Surabaya : Rungkut Menanggal Harapan
Blok O No. 28, Surabaya 60293
Others : Bandung, Makassar, Medan, Balikpapan



(+62) 21 2204 7248
(+62) 31 8786 2593
(+62) 821 8000 1102



www.meditrans.id



info@meditrans.id



[meditrans.id](https://www.facebook.com/meditrans.id)



[meditransid](https://www.instagram.com/meditransid)



Official Partner **MEDITRANS** 

Rivarp Medical Private Limited is a Bangalore based medical devices company, established in 2016. We specialize in R & D, manufacture and sales of high-quality medical devices in therapeutic areas of Interventional cardiology, Spine implants and critical care management.

Rivarp believes in the philosophy of offering best quality medical devices and services to humanity and as a part of this we make sure that each and every product undergoes different quality tests before reaching the operating physician.

Rivarp is accredited by the BSI in (ISO 13485 :2016 & 9001:2015), and passed inspection of Good Manufacturing Practice for Medical Devices. Rivarp Medical's Manufacturing facility is approved for Class D medical devices by the licensing authority Central Drug Standard Control Organization (CDSCO).

Heart Occlusion & Delivery System

STRUCTURAL HEART OCCLUSION SYSTEM / STRUCTURAL HEART DISEASE PRODUCT & DELIVERY SYSTEM.

06

Occlunix® ASD Occluder

12

Occlunix® VSD Occluder

18


Occlunix® PDA Occluder

22


DELIVERY SYSTEM

22

Occlunix® 45/1800 Delivery System

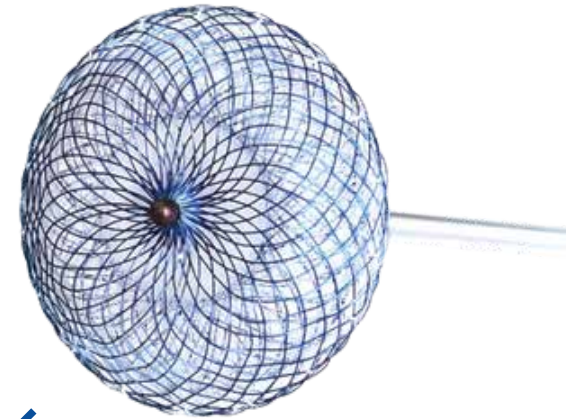


OCCLUNIX® ASD
OCCLUDER

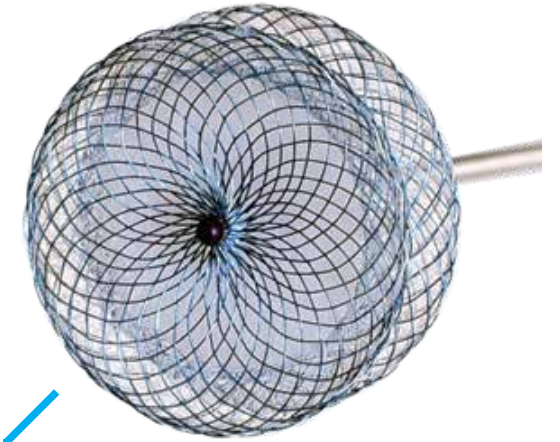

OCCLUNIX® VSD
OCCLUDER



OCCLUNIX® PDA
OCCLUDER

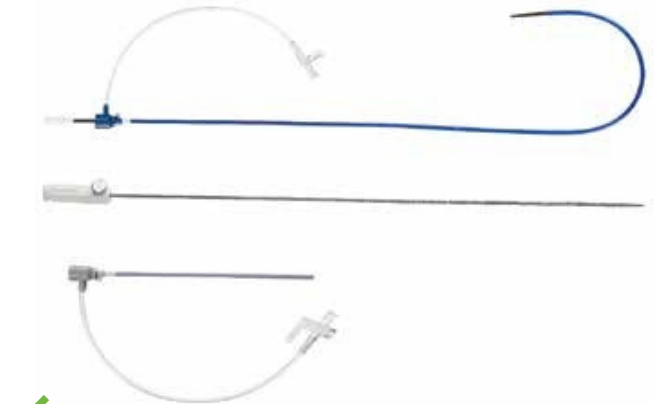


OCCLUNIX® 450 & 1800
DELIVERY SYSTEM



OCCLUNIX® ASD OCCLUDER

SEPTAL OCCLUDER

Occlunix ASD Occluder is a self-expandable, Double-Disc device made from nitinol wire mesh. The frame of the disc provides a foundation for the growth of tissues. The waist that fits the size of defect connects the two discs. Three polyester fabric inserts help in immediate Occlusion.

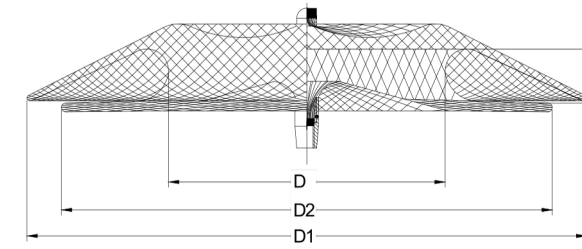

ORDERING INFORMATION

Diameter of discs (mm)

Sl No.	Product Code	Device Size / Waist Dia (D)mm	Left Disc Dia (D1) mm	Waist Length (L)mm	Right Disc Dia (D2)mm	Recommended Sheath Size Occlunix ASD Delivery System (450)
1	RCOA061418	6	18	3.5	14	7F
2	RCOA081620	8	20	3.5	16	8F
3	RCOA101822	10	22	3.5	18	8F
4	RCOA122226	12	26	4	22	9F
5	RCOA142428	14	28	4	24	10F
6	RCOA162630	16	30	4	26	10F
7	RCOA182832	18	32	4	28	10F
8	RCOA203034	20	34	4	30	12F
9	RCOA223236	22	36	4	32	12F
10	RCOA243438	24	38	4	34	12F
11	RCOA263640	26	40	4.5	36	12F
12	RCOA283842	28	42	4.5	38	12F
13	RCOA304044	30	44	4.5	40	14F
14	RCOA324246	32	46	4.5	42	14F
15	RCOA344450	34	50	4.5	44	14F
16	RCOA364652	36	52	4.5	46	14F
17	RCOA384854	38	54	4.5	48	14F
18	RCOA405056	40	56	4.5	50	14F
19	RCOA425258	42	58	4.5	52	14F
20	RCOA445660	44	60	5.5	56	14F
21	RCOA465862	46	62	5.5	58	14F
22	RCOA486064	48	64	5.5	60	14F
23	RCOA506266	50	66	5.5	62	14F




OCCLUNIX® VSD OCCLUDER

SEPTALOCCLUDER

Occlunix VSD Occluder is a self-expandable, double-disc device made from a quality nitinol wire mesh. The two discs are linked together by a cylindrical waist corresponding to the size of the VSD. In order to increase the closing ability, the discs and waist are sewn with three pieces of polyester fabric.

^{technique}
VSD Occluder can be implanted either by venous or arterial catheter

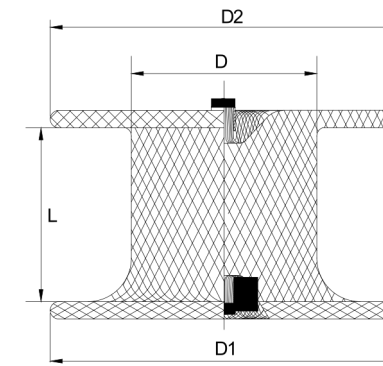

ORDERING INFORMATION

Diameter of discs (mm)


Sl No.	Product Code	Device Size / Waist Dia (D)mm	Left Disc Dia (D1) mm	Right Disc Dia (D2)mm	Waist Length (L)mm	Recommended Sheath Size Occlunix PDA/VSD Delivery System (1800)
1	RCOV050909	5	9	9	3	6 Fr
2	RCOV061010	6	10	10	3	7 Fr
3	RCOV071111	7	11	11	3	7 Fr
4	RCOV081212	8	12	12	3	8 Fr
5	RCOV091313	9	13	13	3	8 Fr
6	RCOV101414	10	14	14	3	8 Fr
7	RCOV121616	12	16	16	3	9 Fr
8	RCOV141919	14	19	19	3	10 Fr
9	RCOV162121	16	21	21	3	10 Fr
10	RCOV182323	18	23	23	3	12 Fr
11	RCOV051111	5	11	11	7	6 Fr
12	RCOV061212	6	12	12	7	7 Fr
13	RCOV071313	7	13	13	7	7 Fr
14	RCOV081414	8	14	14	7	8 Fr
15	RCOV091515	9	15	15	7	8 Fr
16	RCOV101616	10	16	16	7	8 Fr
17	RCOV121818	12	18	18	7	9 Fr
18	RCOV142020	14	20	20	7	10 Fr
19	RCOV162222	16	22	22	7	10 Fr
20	RCOV182424	18	24	24	7	12 Fr



OCCLUNIX® PDA OCCLUDER

Occlunix PDA Occluder is a self-expandable device made from a quality nitinol wire mesh in a T-cone shaped design. The design allows the device to conform eventually to the ducts shape for optimal occlusion.

The pieces of polyester fabric sewn inside to help immediate occlusion and faster reendothelialization

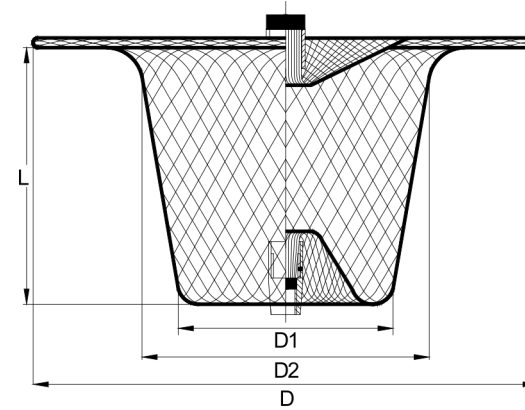


ORDERING INFORMATION

Diameter of discs (mm)

Sl No.	Product Code	Retention Skirt (D) mm	Proximal Diameter (D1) mm at screw	Distal Diameter (D2)mm waist diameter	Device Length (L)mm	Recommended Sheath Size Occlunix PDA/VSD Delivery System (1800)
1	RCOP040610	10	4	6	7	6F
2	RCOP060812	12	6	8	7	6F
3	RCOP081014	14	8	10	7	7F
4	RCOP101216	16	10	12	7	8F
5	RCOP121420	20	12	14	7	8F
6	RCOP141622	22	14	16	8	8F
7	RCOP161824	24	16	18	8	9F
8	RCOP182026	26	18	20	9	10F
9	RCOP202228	28	20	22	9	12F
10	RCOP222430	30	22	24	10	12F

OCCLUNIX® 450 & 1800 DELIVERY SYSTEM

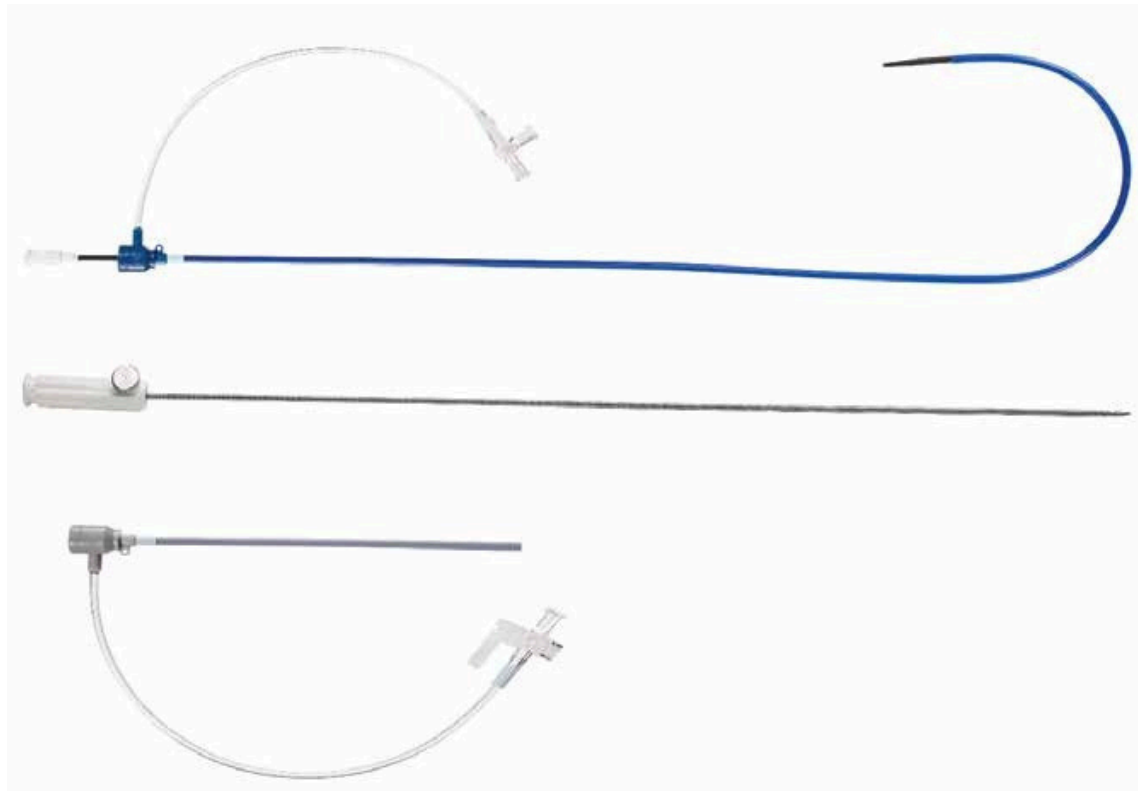

Occlunix 450 Delivery Sheath ASD

Sl No.	Product Code	Sheath Size (Fr)	Length (mm)
1	RCODS0745	7F	856
2	RCODS0845	8F	856
3	RCODS0945	9F	856
4	RCODS1045	10F	856
5	RCODS1245	12F	856
6	RCODS1445	14F	856

Occlunix 450 Delivery System consists of a Delivery Sheath, Dilator, Loader, Hemo Static Valve and delivery cable for the delivery of ASD Occluder.

Occlunix 1800 Delivery Sheath VSD & PDA

Sl No.	Product Code	ID (Delivery Sheath)	length (mm)
1	RCODS06180	7F	856
2	RCODS07180	8F	856
3	RCODS08180	9F	856
4	RCODS09180	10F	856
5	RCODS10180	12F	856
6	RCODS12180	14F	856

Occlunix 1800 Delivery System Consists of a Delivery Sheath, Dilator, Loader, Hemo Static Valve and delivery cable for the delivery of PDA & VSD Occluder.